



Mobile PERS – On The Go

On the Go

Designed to provide peace of mind, On the Go mobile is versatile, small, and easy to use delivering 24/7 fast access to help at the press of a button – both in and outside the home




Key features:

- 4G LTE Cat1 device delivers a faster wireless communication speed and keeps users connected without the need of a cellular contract or landline service
- GPS and Wi-Fi location technology enables us to send assistance to the users' exact location
- Built-in speaker allows users to speak directly to a Trained Care Specialist when assistance is needed
- Long-lasting battery life reduces the need for daily charging
- Waterproof pendant can be safely worn in the shower or tub
- Optional fall detection automatically connects to the response center, even if the user cannot press the help button

What's Included In The Box


- Mobile Device
- Charging Cradle
- Belt Clip
- Neck Lanyard
- Quick start guide

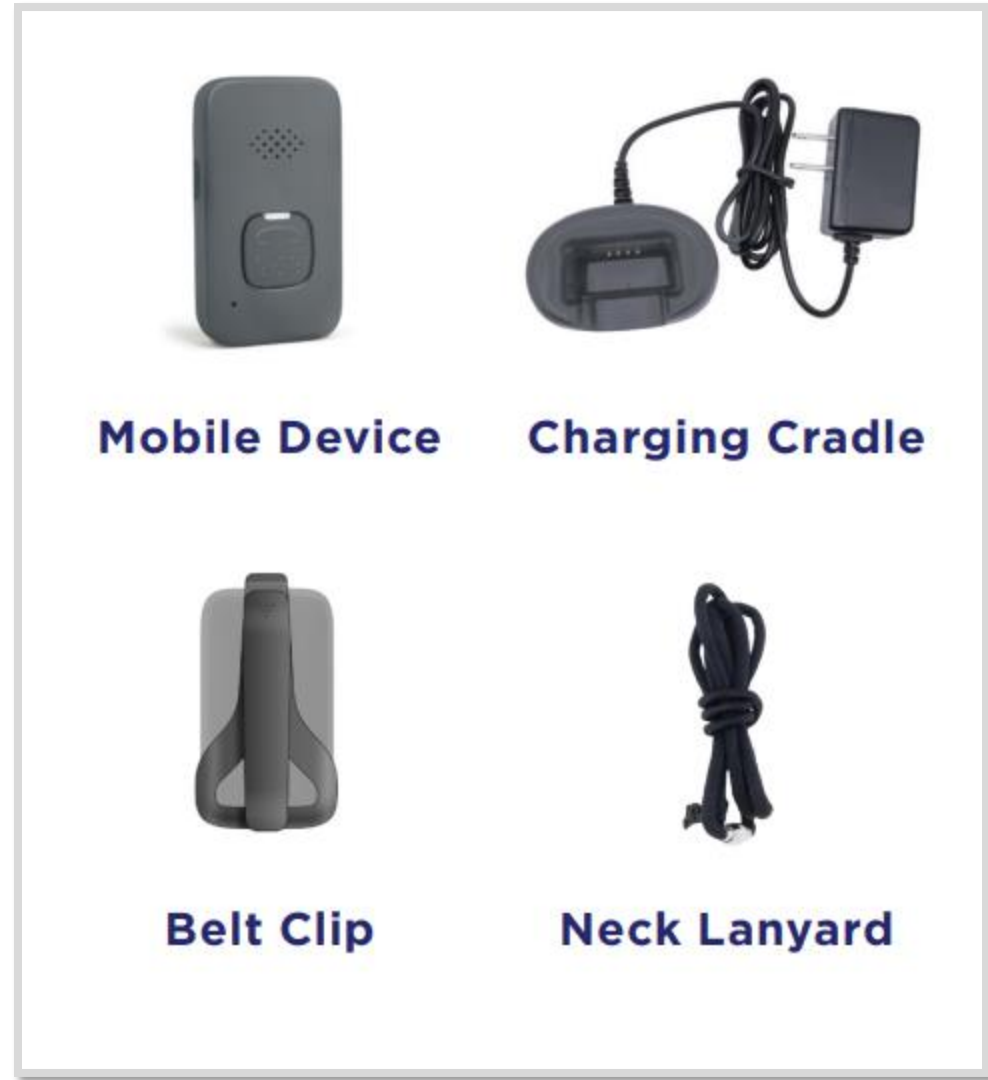

Device Dimensions





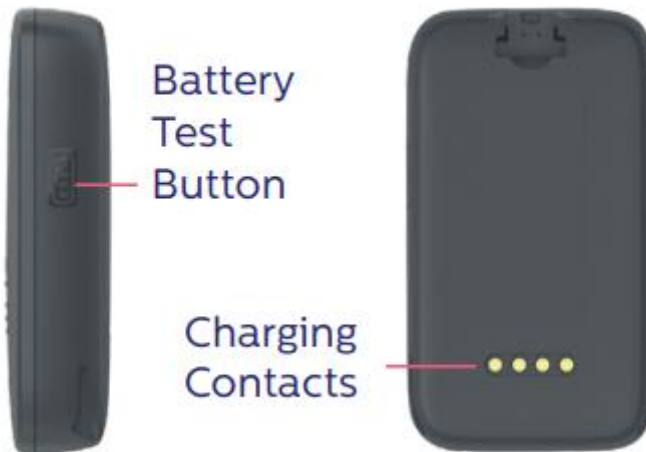


System Overview



Status Light Chart

Battery Status

- Confirm battery status, press the battery test button
- Device will announce “*Battery Low*” or “*Battery Okay*”
- Device is designed to be on at all times

<u>Status Light Chart</u>		
Status Light	Definition	Visual
Solid Red	Fully Charged	
Flashing Red	Charging	
Solid Blue	Placing a Call	
Flashing Blue	Normal Operation	
No Light On	Device Battery is Depleted	

Mobile PERS – On The Go

Tips and Reminders

- The battery life can last 4 to 6 days
- Performance of GPS can be affected by a wide range of factors
- Device is waterproof but should not be submerged for extended periods
- Lanyards are designed to breakaway under certain conditions
- It is essential to test monthly to ensure it is functioning properly



Lifeline